

100% Ethylene Oxide sterilization

ANDERSEN PRODUCTS  
www.sterility.com

NOW  
YOU  
CAN

## The Scope Sterilization Solution

*The most popular low temperature sterilization system in North America keeps getting better! Using less than 20 grams of Ethylene Oxide per cycle, Anprolene® gently sterilizes your delicate scope that heat and moisture can damage. Our system's unique bag and ampoule method allows for load size flexibility. Easy installation requires only 110V outlet and a dedicated exhaust line to the outside.*

### ► Unmatched Sterilization

**Flexibility:** Anprolene's 100% EtO System sterilizes at room temperature without a vacuum or water injection. Protect your most delicate instruments, scopes, drills, fiber optics, cameras, and plastic or rubber items... Literally everything but food, drugs or liquids can be sterilized with Anprolene.

### ► Sterilization and Aeration in the same cabinet:

The Anprolene system utilizes an advanced purge system that aerates products for two hours at the end of the cycle. Active aeration continues indefinitely until products are removed.

### ► Unique Load Size Flexibility:

The Anprolene system uses unique plastic sterilization liner bags and individual gas ampoules. We offer two different ampoule/bag combinations to match the size of your sterilization load. You never use more gas than you have to!

### ► Proven Reliability:

Anprolene has been the choice of healthcare professionals for over 40 years. Anprolene is used in over 5,000 facilities worldwide.



**UPGRADE**  
to bags made for the  
EDSS for the same price!  
Part: BVTAN74.87BV

Details on back

### AN74ix Starter kit

~~RETAIL~~ SPECIAL  
~~\$4060.40~~ \$3999.00

#### STARTER KIT INCLUDES

Units	Item	Reorder	Quantity	Retail
25 each	5 cc ampoules and sterilization liner bags	AN71	1 case	\$145.00
20 each	20 cc ampoules and sterilization liner bags	AN79	1 case	\$220.00
50 each	Humidichip relative humidity stabilizer chips	AN1071	2 jars	\$ 17.00
1 each	Humiditube holder	AN1072	1 each	\$ 1.00
50 each	Dosimeter chemical indicators	AN87	2 boxes	\$ 57.40
400 each	EtO exposure indicator strips - chemical indicators	AN85	2 boxes	\$ 55.00
22 each	Steritest biological indicators	AN80	2 boxes	\$100.00
2 rolls	3" Seal & Peel (200 feet per roll)	AN830	1 box	\$ 70.00
1 roll	5" Seal & Peel (200 feet per roll)	AN850	1 box	\$ 55.00
1 roll	7" Seal & Peel (200 feet per roll)	AN870	1 box	\$ 75.00
1 each	Impulse sealer	AN90	1 each	\$130.00
1 each	AirScan test kit (STEL and TWA)	AN93	1 kit	\$ 85.00

BIO  
VISION Technologies, LLC  
RECOMMENDED

**ANPRO** ANDERSEN  
PRODUCTS

the sterility experts

www.sterility.com  
Andersen Products, Inc.  
Health Science Park • 3202 Caroline Drive  
Haw River, NC 27258-9564 USA  
800-523-1276 • fax 336-376-8153

www.andcal.com  
H. W. Andersen Products of California, Inc.  
Health Science Park • 3151 Caroline Drive  
Haw River, NC 27258-9789 USA  
800-524-3455 • fax 336-376-3088

Distributed in Europe by:  
H.W. Andersen Products, Ltd.  
Davy Road • Clacton-on-Sea  
Essex CO15 4XA UK  
1255-428-328


Manufactured by:  
Andersen Sterilizers, Inc.  
3154 Caroline Drive  
Haw River, NC 27258 USA

## Sterilizer technical specifications

- ▶ **Longer Cabinet Capacity:** The AN74ix has an extended cabinet to allow for longer items and larger loads. The interior cabinet is 29.25" in length.
- ▶ **Ease of Installation:** New High volume and Low volume pressure pumps simplify installation with only a 7/8" Reinforced PVC hose connection to the outside. The ventilation and control panel are built in on top of the sterilizer cabinet.
- ▶ **Sterility Assurance:** Convenient tests are available for each cycle to assure success. Steritest® is a Biological Indicator (BI), the classic test of sterility. BI's require a culture period. Immediate results are provided by Dosimeter®, a Chemical Indicator, and a rigorous test for every cycle. Anprolene Exposure Indicators are used for each item within a "sterilization load" to ensure individual exposure.

### Specifications: AN74ix

Model AN74.87  
Interior dimensions: 29.25"(L) x 13.75"(W) x 13.75"(H)  
Exterior dimensions: 30.50"(L) x 14" (W) x 20.5"(H)  
Total weight: 60.5 lbs.  
Power Supply: 115v +/- 10%, 60 Hz/150 Watts  
Installation Category II  
Manufactured by: Andersen Sterilizers, Inc.  
Distributed by: Andersen Products, Inc.

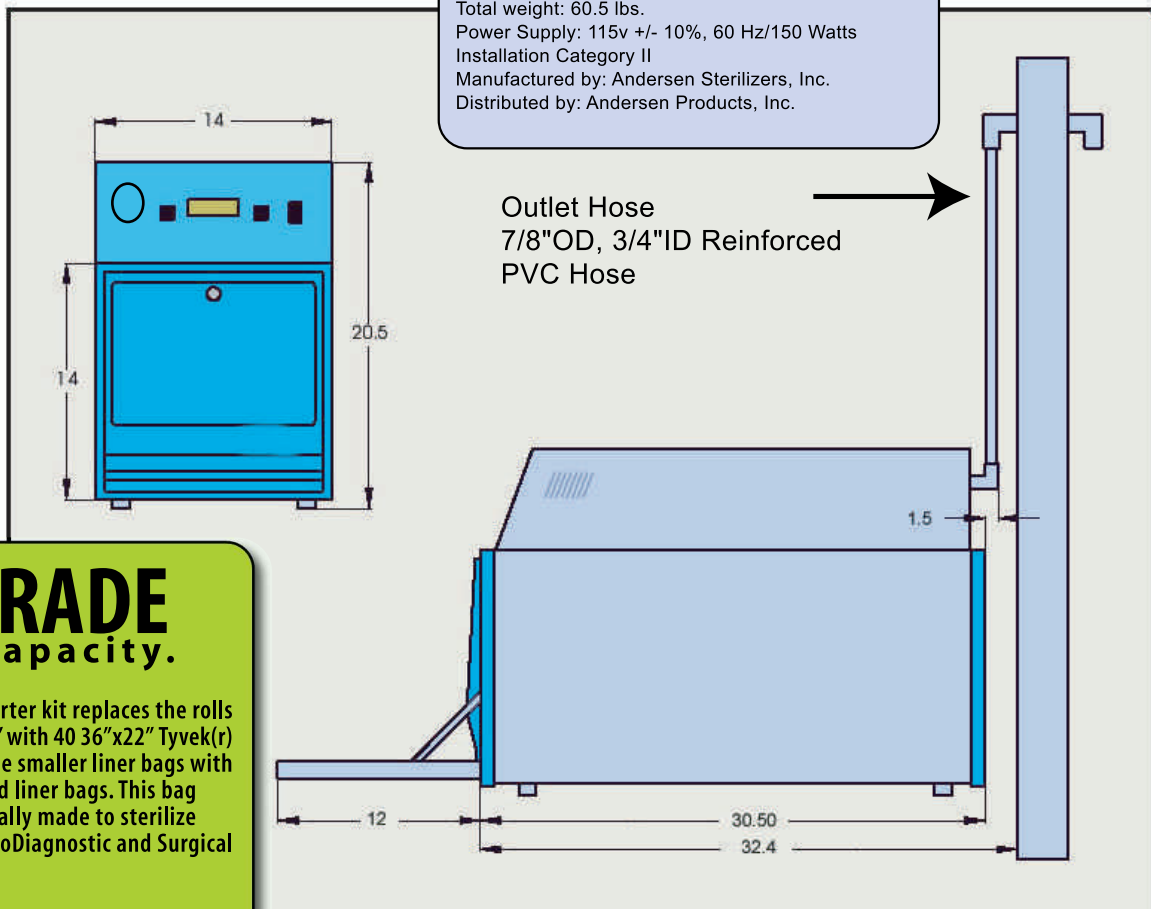

## UPGRADE your capacity.

This special starter kit replaces the rolls of "Seal & Peel" with 40 36"x22" Tyvek(r) pouches and the smaller liner bags with 10 22"x42" load liner bags. This bag system is specially made to sterilize BioVision's EndoDiagnostic and Surgical Suite (EDSS).

Andersen Part: BVTAN74.87BV

Steritest® is a registered trademarks of Andersen Products, Inc.

Call for further information about Anprolene® Sterilization Kits and prices.  
Call 800-523-1276