

BioVision Technologies


Endoblade

4 IN 1 ENDOSURGICAL INSTRUMENT

The BioVision Veterinary Endoscopy EndoBlade features:


- ✓ Reduced instrumentation count
- ✓ Speeds procedure time
- ✓ Helps maintain insufflation
- ✓ Provides excellent hemostasis
- ✓ Compatible with the BioVision SurgAssist

The 4 Functions of the EndoBlade

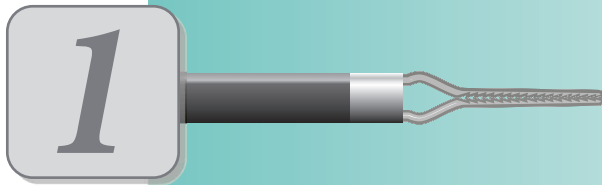


Probe

The blunt tip allows for atraumatic manipulation of Tissues.

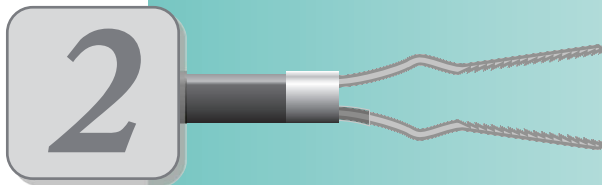

Grasping Forceps

Atraumatic teeth provide a firm grasp.

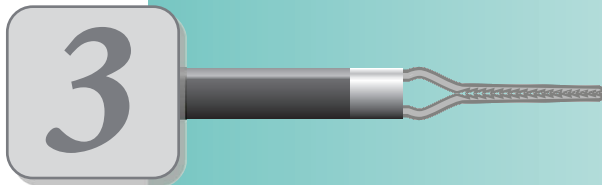

Bipolar Cautery

Compression cautery is accomplished by passing RF energy between the forcep paddles.



Bipolar Cutting with Cautery

Cutting with cautery is accomplished by passing RF energy between the blade and forcep paddles.

Specifications:

O.D. - 5.16mm

W.L. - 430mm

Generator connection - dual 4mm "banana" plugs

Part No.: BVTINST4-00001



BioVision Veterinary Endoscopy, 221 Corporate Circle Suite H, Golden, Colorado 80401

biovisiontech.com

BVT93065-00001-A
120904